

Cable Wiring Guide

1263B - TECSOURCE Cable, 10A, 2m, Pigtailed

Length: 2 meter


Instrument End

Connector: DB-15 Male

Label: Line 1: 1263B
Line 2: TECSOURCE, 10A, 2m

Pin	Pair	Function / Description	Gauge
1	RED	TE+	20 AWG
2	RED/WHT	TE+	20 AWG
3	BLK	TE-	20 AWG
4	BLK/WHT	TE-	20 AWG
5	Shield	Earth Ground	22 AWG
6	Shield	Earth Ground	22 AWG
7	GRN	Sensor+	24 AWG
8	WHT	Sensor-	24 AWG
9	RED/YEL	TE+	20 AWG
10	BLK/YEL	TE-	20 AWG
11	YEL	Fan+	24 AWG
12	BLU	Fan-	24 AWG
13	Short to pin 12	(ID Pin, Male Side Only)	Short to pin 12
14	ORG	R Sens+	24 AWG
15	BRN	R Sens-	24 AWG
Shell	Shield	Earth Ground	22 AWG

DB-15 Male

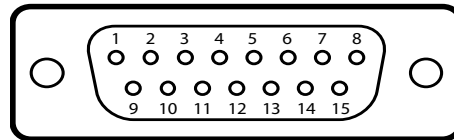

Device/Mount End

Connector: None, bare wire (pigtail is 3")